

SuperCross[™]
microcatheters

0.014" STRAIGHT TIP & FT



SUPERIOR CROSSABILITY

with exceptional deliverability

SuperCross™ microcatheters



SuperCross microcatheters offer the right flexibility for superior deliverability and guidewire support.

STAINLESS STEEL BRAID

for flexibility, pushability and kink resistance


FULL-LENGTH LUMEN

facilitates subselective delivery of contrast and includes a PTFE inner layer for easy guidewire exchange

HYDROPHILIC COATING

for smooth delivery through tortuous vessels

FLEXIBLE TIP (FT) VERSION PROVIDES EXCEPTIONAL DELIVERABILITY

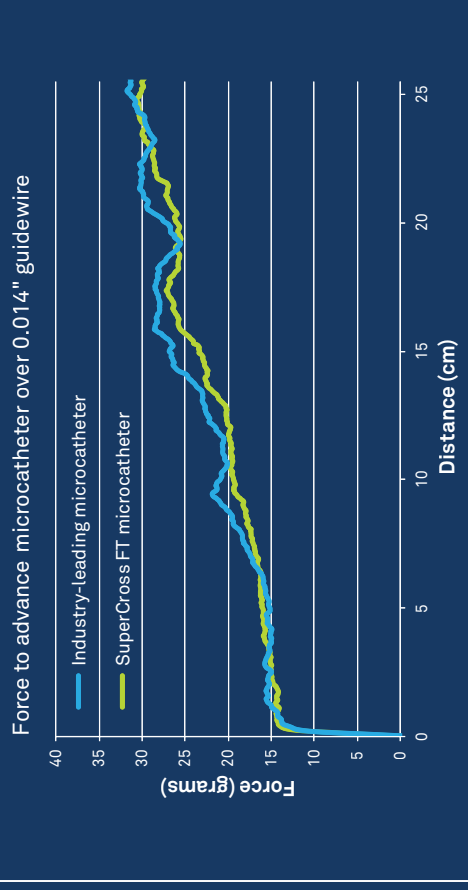

STRAIGHT TIP VERSION

for superior guidewire tracking

FLEXIBLE TIP (FT) VERSION

for exceptional deliverability

Embedded platinum/iridium marker band

creates brilliant radiopacity for confirming device placement




ORDERING INFORMATION

The SuperCross microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

Model	Type	Guidewire Compatibility	Distal Tip	Working Length	Shaft Construction	Distal O.D.	Proximal O.D.	Distal I.D.	Proximal I.D.	Hydrophilic Coating
5300	OTW	0.014"	Straight	130cm	Braided	1.8F (0.61mm)	2.5F (0.84mm)	0.017" (0.43mm)	0.021" (0.53mm)	Distal 40cm
5301	OTW	0.014"	Straight	150cm	Braided					
5340	OTW	0.014"	FT	130cm	Braided	1.8F (0.61mm)	2.5F (0.84mm)	0.017" (0.43mm)	0.021" (0.53mm)	Distal 40cm
5341	OTW	0.014"	FT	150cm	Braided					

Packaged in quantities of 1 unit per box.
Flexible Tip = FT

STRAIGHT - FT

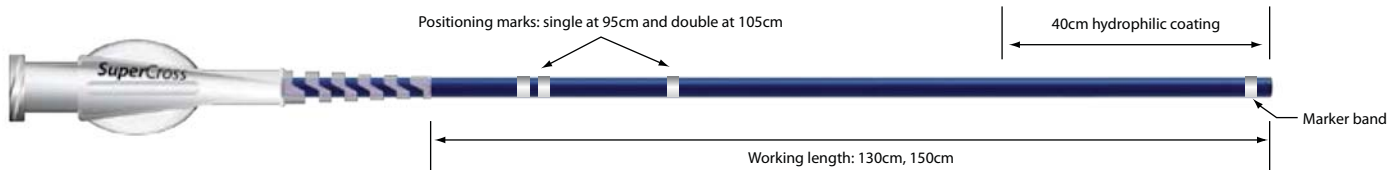

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

SuperCross is a trademark of Vascular Solutions, Inc.

©2014 Vascular Solutions, Inc. All rights reserved. ML2940 Rev. A 11/14



Vascular Solutions, Inc.
6464 Sycamore Court North
Minneapolis, Minnesota 55369 USA

Customer Service:
United States: 888.240.6001
International: (001) 763.656.4298
customerservice@vasc.com
www.vasc.com

0.014" AT

SuperCross[™]
microcatheters




SUPERIOR CROSSABILITY

in tortuous anatomy



SuperCross microcatheters provide support in tortuous anatomy and navigating bifurcated vessels.

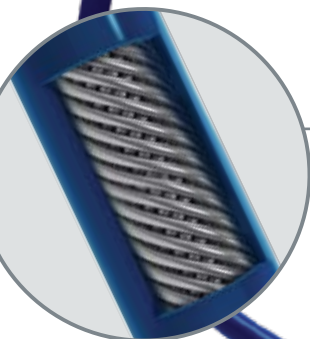

DUAL COIL DESIGN

SuperCross Angled Tip OTW provides excellent torque response, flexibility, pushability and kink resistance

FULL-LENGTH LUMEN

facilitates subselective delivery of contrast


HYDROPHILIC COATING

on the distal 80cm



THE RIGHT TORQUE RESPONSE FOR TORTUOSITY

The new dual coil design of the SuperCross Angled Tip OTW significantly improves torque response

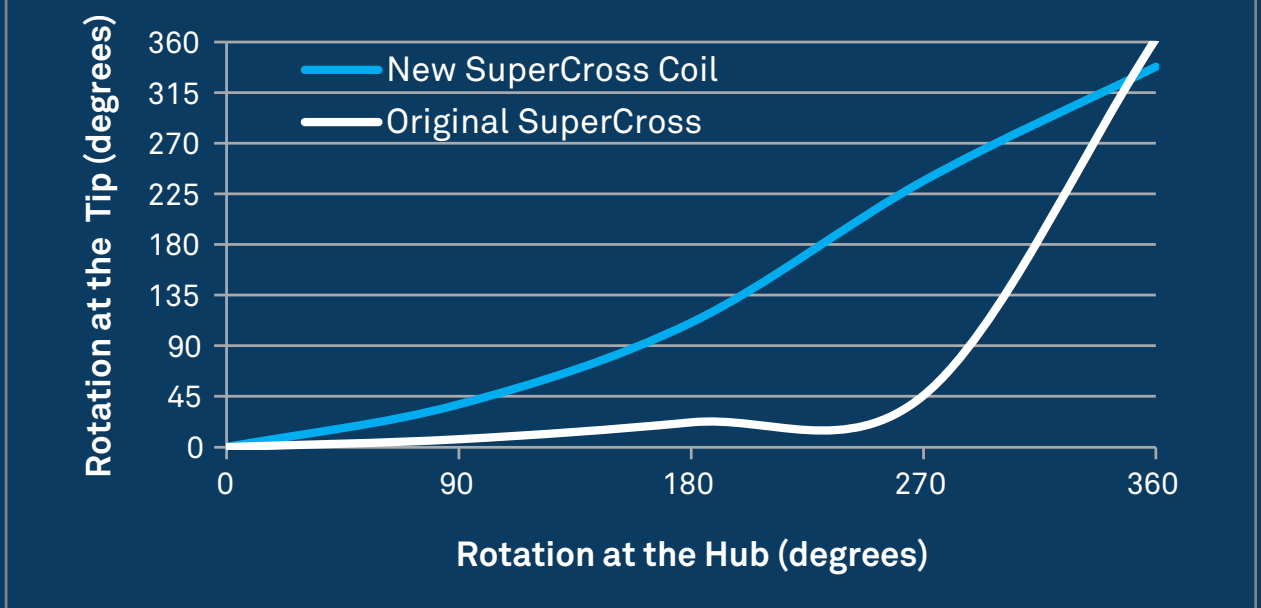


Testing completed by Vascular Solutions, Inc. Data on file.

MULTIPLE ANGLE OPTIONS

Embedded platinum/tungsten coil provides enhanced visibility along the entire angled tip

XT= 90° extended tip for secure cannulation

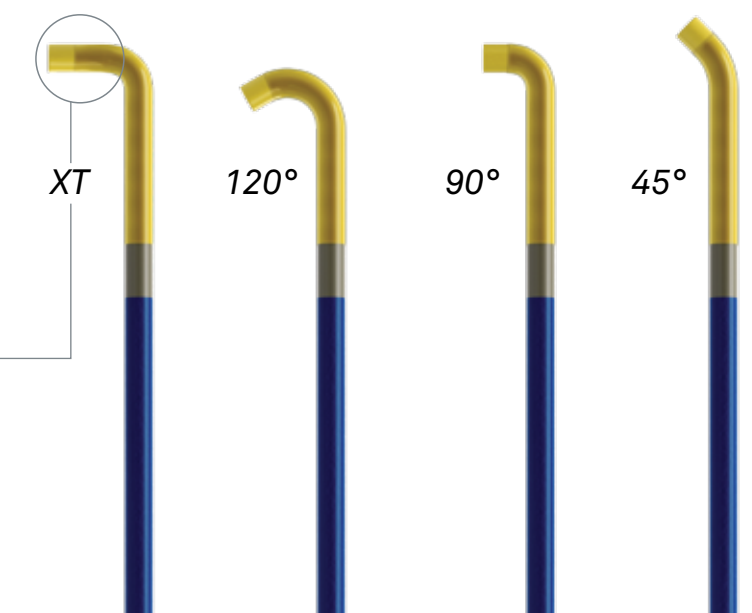


ORDERING INFORMATION

The SuperCross microcatheter is intended to be used in conjunction with steerable guidewires to access discrete regions of the coronary and/or peripheral vasculature. It may be used to facilitate placement and exchange of guidewires and other interventional devices and to subselectively infuse/deliver diagnostic and therapeutic agents.

Model	Type	Guidewire Compatibility	Distal Tip	Working Length	Shaft Construction	Distal O.D.	Proximal O.D.	Distal I.D.	Proximal I.D.	Hydrophilic Coating
5302	OTW	0.014"	45°	130cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5303	OTW	0.014"	45°	150cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5304	OTW	0.014"	90°	130cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5305	OTW	0.014"	90°	150cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5306	OTW	0.014"	120°	130cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5307	OTW	0.014"	120°	150cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5308	OTW	0.014"	90° XT	130cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm
5309	OTW	0.014"	90° XT	150cm	Coiled	2.4F 0.79mm / 0.031"	3.2F 1.07mm / 0.042"	0.017" (0.43mm)	0.018" (0.46mm)	Distal 80cm

Packaged in quantities of 1 unit per box.

ANGLED TIP

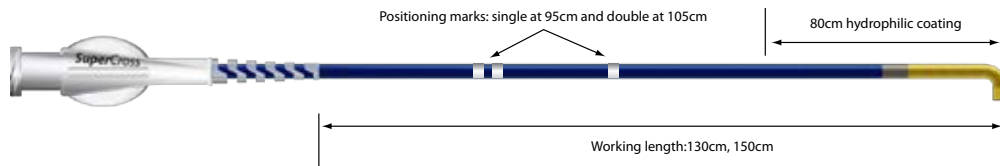


ALSO AVAILABLE IN STRAIGHT - FLEXIBLE TIP

Model	Type	Guidewire Compatibility	Distal Tip	Working Length	Shaft Construction	Distal O.D.	Proximal O.D.	Distal I.D.	Proximal I.D.	Hydrophilic Coating
5300	OTW	0.014"	Straight	130cm	Braided	1.8F 0.61mm / 0.024"	2.5F 0.84mm / 0.033"	0.017" (0.43mm)	0.021" (0.53mm)	Distal 40cm
5301	OTW	0.014"	Straight	150cm	Braided	1.8F 0.61mm / 0.024"	2.5F 0.84mm / 0.033"	0.017" (0.43mm)	0.021" (0.53mm)	Distal 40cm
5340	OTW	0.014"	FT	130cm	Braided	1.8F 0.61mm / 0.024"	2.5F 0.84mm / 0.033"	0.017" (0.43mm)	0.021" (0.53mm)	Distal 40cm
5341	OTW	0.014"	FT	150cm	Braided	1.8F 0.61mm / 0.024"	2.5F 0.84mm / 0.033"	0.017" (0.43mm)	0.021" (0.53mm)	Distal 40cm

Packaged in quantities of 1 unit per box.
Flexible Tip = FT

STRAIGHT - FT

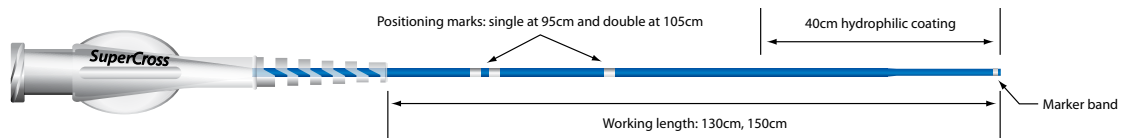

Please see the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.

CAUTION: Federal Law (U.S.A.) restricts this device to sale by or on the order of a physician.

SuperCross is a trademark of Vascular Solutions, Inc.

©2016 Vascular Solutions, Inc. All rights reserved. ML2938 Rev. E 07/16



Vascular Solutions, Inc.
6464 Sycamore Court North
Minneapolis, Minnesota 55369 USA

Customer Service:
United States: 888.240.6001
International: (001) 763.656.4298
customerservice@vasc.com
www.vasc.com